

DATA SHEET

eFISH NMYC/ 2q11Dual ColorProbe

Catalog No.**FP043-10XE-100µl - 10 test****FP043-20XE-200µl - 20 test**

Doc No: 932-FP043E Rev: D

Date of Release: 05-Aug-2020

Material Provided: One vial of eFISH probe in hybridization buffer (RTU).

Recommended detection system (Not supplied):

Either of the following detection system is recommended depending on the automation/manual platform used:

eFISH Kit	Cat #	Description
eFISH Histo	DF-500-20XE	Automation
eFISH Cyto	DF-510-20XE	Automation

Intended Use:

The BioGenex eFISH NMYC/2q11 Dual Color probe is designed to be used for the detection of the human NMYC gene as well as chromosome 2 specific sequences in formalin-fixed, paraffin-embedded tissue or cells by fluorescence in situ hybridization (FISH).

BioGenex eFISH NMYC/2q11 Dual Color probe comes in hybridization buffer. The probe contains green-labeled polynucleotides (Green: excitation at 503 nm and emission at 528 nm, similar to FITC), which target the NMYC gene, and orange-labeled polynucleotides (Orange: excitation at 547 nm and emission at 572 nm, similar to rhodamine), which target sequences of chromosome 2 in the chromosomal region 2q11

Summary and Explanation

Fluorescence *in situ* hybridization (FISH) is a robust technique of cytogenetic used for the detection of chromosomal aberrations, presence or absence of specific DNA sequence in native context. In this technique fluorescent probes bind to the target sequence of DNA in chromosome. High specificity and sensitivity coupled rapid and an accurate result has proven role of FISH in both research and diagnosis of solid tumor and hematological malignancies. As technique of cancer cytogenetics, FISH, can be used to identify genetic aberrations viz., deletions, amplification and translocation in tissue sections or within individual cells. FISH is also used for use in genetic counseling, medicine, and species identification. FISH can also be used to detect and localize specific RNA targets in cells,

circulating tumor cells, and tissue samples^{1,2,3,4,5}.

In FISH procedure, fixed tissue sections are pretreated to expose target DNA or mRNA sequences. An appropriately labeled probe is hybridized to the exposed target DNA or mRNA sequences in the cells. Subsequent stringent washing steps remove any probe that is non-specifically bound to the tissue section. Subsequently slides are mounted using DAPI/antifade and can be visualized under fluorescence microscope using appropriate filter set.

Principles of the Procedure

In Situ hybridization (ISH) allows the detection and localization of definitive nucleic acid sequences directly within a cell or tissue. High specificity is ensured through the action of annealing of fluorescence probe nucleic acid sequence to complementary target nucleic acid sequence. ISH techniques can be used to identify genetic aberrations like deletions, amplification, and translocation in tissue sections or within individual cells.

Storage and Handling

The BioGenex eFISH NMYC/2q11 Dual Color Probe must be stored at 2-8°C protected from light and is stable through the expiry date printed on the label.

Specimen Collection and Slide Preparation

Tissues fixed in 10% (v/v) formalin are suitable for use prior to paraffin embedding and sectioning.

FISH Staining procedure

- (a) The BioGenex eFISH probes are supplied in hybridization buffer and used without further dilution.
- (b) Protocol:

Please refer to the eFISH probe specific instruction/protocol for automated or semi-automated FISH processing platform (Xmatrx[®]-Infinity, Xmatrx[®]-Nano and Xmatrx[®]-mini.

Further processing, such as washing and counter-staining, can be completed according to the user's needs. For a particularly user-friendly performance, we recommend the use of a BioGenex eFISH kit.

These systems were also used for the confirmation of appropriateness of the BioGenex eFISH NMYC/2q11 Dual Color Probe.

Disclaimer: The above information is provided for reference only. Each end-user is responsible for developing and validating optimal testing conditions for use with this product.

Troubleshooting

Contact BioGenex Technical Service Department at **1-800-421-4149** or your local distributor to report unusual staining.

Expected Results

The BioGenex eFISH NMYC/2q11 Dual Color Probe is a mixture of a green fluorochrome direct labeled NMYC probe hybridizing to the human NMYC gene in the chromosomal region 2p24.3 and an orange fluorochrome direct labeled 2q11 probe specific for the AFF3 (AF4/FMR2 family, member 3) gene region in 2q11.2. Due to cross-hybridizations of chromosome 2 alpha satellites to other centromeric regions, probes specific for 2q11 are frequently used for chromosome 2 copy number detection.

In a normal interphase nucleus, two orange and two green signals are expected. In a cell with amplification of the MYC gene locus, multiple copies of the green signal or green signal clusters will be observed.

However, we recommend the use of a control sample in which the 2q11 status is known to judge the specificity of the signals with each hybridization reaction.

Care should be taken not to evaluate overlapping cells, in order to avoid false results, e.g. an amplification of genes. Due to decondensed chromatin, single FISH signals can appear as small signal clusters. Thus, two or three signals of the same size, separated by a distance equal to or less than the diameter of one signal, should be counted as one signal.

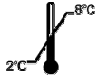






Limitations of the Procedure

Correct treatment of tissues prior to and during fixation, embedding, and sectioning is important for obtaining optimal results. Inconsistent results may be due to variations in tissue processing, as well as inherent variations in tissue. The results from *in situ* hybridization must be correlated with other laboratory findings.

Bibliography

1. Gall, J. G. and Pardue, M. L. (1969). *Proc. Natl. Acad. Sci. USA* 63,378 -383.
2. Rudkin, G. T. and Stollar, B. D. (1977). *Nature* 265,472-473.
3. Hougaard, D. M., Hansen, H. and Larsson, L. I. (1997). *Histochem. Cell Biol.* 108,335 -344.
4. Bauman, J. G., Wiegant, J., Borst, P. and van Duijn, P. (1980). *Exp. Cell Res.* 128,485 -490.
5. O'Connor, C. (2008) *Nature Education* 1(1):171.

6. David N. Shapiro, Marcus B. Valentine, Susan T. Rowe, Anne E. Sinclair, Jack E. Sublett, W. Mark Roberts, and A. Thomas Look. (1993). Detection of N-myc Gene Amplification by Fluorescence in Situ Hybridization Diagnostic Utility for Neuroblastoma. *American journal of Pathology*, Vol. 142, No. 5.
7. Miao Wang¹, Chunju Zhou², Rongqin Cai¹, Yong Li¹ and Liping Gong. (2013). *Diagnostic Pathology*, 8:5
8. Prasad Mathew, Marcus B Valentine, and A Thomas Look. (2001). *Vol 3, No 2*, pp105-109.

	Temperature Limitation		In Vitro Diagnostic Medical Device
	Use By Date		Batch Code
	Non-Sterile		Consult Instructions for Use
	Representative in the European Community	