

Anti-Cytomegalovirus (CMV) [DDG9/CCH2]

Catalog No.	Description
AM997-5M	6 ml of Ready-to-Use Antibody for use with BioGenex Super Sensitive™ Detection Systems OR equivalent detection system
AM997-10M	10 ml of Ready-to-Use Antibody in a barcode labeled vial for use with BioGenex Super Sensitive™ Detection Systems and i6000™ Automated Staining Systems
MU997-UC	1 ml of Concentrated Antibody for use with BioGenex Super Sensitive™ Detection Systems OR equivalent detection system
MU997-5UC	0.5 ml of Concentrated Antibody for use with BioGenex Super Sensitive™ Detection Systems OR equivalent detection system
AX997-YCD	Ready-to-Use Antibody in Barcode labeled vial for use on the Xmatrx® Elite/Ultra Staining System, 160 tests
AX997-50D	Ready-to-Use Antibody in Barcode labeled vial for use on the Xmatrx® Elite/Ultra Staining System, 50 tests

Clone	Species	Ig Class
DDG9/CCH2	Mouse	IgG2a (DDG9), IgG1 (CCH2), kappa

Intended Use

For Research Use. This antibody is designed for the specific localization of Cytomegalovirus (CMV) in formalin-fixed, paraffin-embedded (FFPE) tissue sections. Evaluation must be performed by a qualified pathologist.

Summary and Explanation

Cytomegalovirus (CMV) is a member of the family Herpesviridae, which includes Epstein-Barr virus, herpes simplex virus types 1 and 2, varicellazoster virus, and human herpes virus 6. CMV is the most common identified cause of congenital infection. Found in several body fluids including saliva, urine, breast milk, cervical secretions, blood, and semen, CMV can be transmitted in a variety of ways, such as blood transfusion and organ transplantation. The typical morphology of infected tissue includes enlarged cells (often to a diameter of 40 μm) with intranuclear, and on occasion, cytoplasmic inclusions. The intranuclear inclusions may measure 17 μm in diameter and are usually set off from the nuclear membrane by a clear halo.

Storage and Handling

Store at 2-8°C. Fresh dilutions, if required, should be prepared prior to use and are stable and steady for up to one day at room temperature (20-26°C). Diluted antibody preparations can be refrigerated or frozen for extended shelf life.

Principles of the Procedure

Antigen detection by immunohistochemistry (IHC) is a two-step process wherein the primary antibody binds to the antigen of interest and that binding is detected by a chromogen. The [primary antibody](#) may be used in IHC using manual techniques or BioGenex Automated Staining System. Positive and negative controls should always be run simultaneously with all patient specimens.

Reagents Provided

Mouse Monoclonal Antibody Cytomegalovirus (CMV) is affinity purified and diluted in PBS, pH 7.2, containing 1% BSA and 0.09% sodium azide

Dilution of Primary Antibody

BioGenex Ready-to-Use antibodies have been optimized for use with the recommended BioGenex Detection System and should not require further dilution.

BioGenex concentrated antibodies must be diluted in accordance with the recommended protocol when used with the recommended BioGenex Detection System.

Recommended Protocol

Refer to the following table for conditions specifically recommended for this antibody. Refer to the BioGenex website for guidance on specific staining protocols or other requirements.

Parameter	BioGenex Recommendations
Control Tissue	CMV infected lung tissue as available with Biogenex FB-997M* & FG-997M*
Recommended Dilution for Concentrated Antibody	1:10-25 in HK941
Recommended Pretreatment (Manual/i6000)**	EZ-AR2 (HK522-XAK)
Recommended Pretreatment (Xmatrx)	EZ-AR2 Elegance (HX032-YCD)
Antibody Incubation (Manual/i6000)	30-60 Min at RT
Antibody Incubation (Xmatrx)	45-60 Min at 25°C
Detection System for Manual, Xmatrx & i6000 systems***	Use BioGenex Two-Step OR One-Step Super Sensitive™ Polymer-HRP IHC Detection System/DAB; see p. 2 for more information

*FB: positive control micro chamber slides, FG: positive control microscopic slides. Xmatrx requires micro chamber slides.

**Pretreatment times will vary based on individual microwave power.

Category	Antibodies	Revision No.	D
Document No.	932-997M-EN	Release Date	22-Nov-2021

***For automation systems (Xmatrix-Elite, Xmatrix-Ultra & i6000 Diagnostics), refer to the factory protocols provided with the instrument.

Detection System	Two-Step HRP Kit	One-Step HRP Kit	Link and Label Kit
Manual	QD440-XAKE (1000 Test)	QD630-XAKE (1000 Test)	QP300-XAKE (1000 Test)
	QD430-XAKE (1000 Test)		
	QD420-YIKE (500 Test)	QD620-XAKE (500 Test)	QP900-9LE (500 Test)
	QD400-60KE (60 Test)		
Xmatrix - Automation	QD550-YCDE (200 Test)	QD610-YADE (200 Test)	N/A
i6000 - Automation	QD410-YAXE (200 Test)	QD610-YAXE (200 Test)	N/A
For more information, visit www.biogenex.com .			

Precautions

This product contains sodium azide at concentrations of less than 0.1%. Sodium azide is not classified as a hazardous chemical at the product concentrations, but proper handling protocols should be observed. For more information, a Safety Data Sheet (SDS) for sodium azide is available upon request. Dispose of unused reagents according to Local, State and Federal Regulations. Wear suitable Personal Protective Equipment, do not pipette reagents by mouth, and avoid contact of reagents and specimens with skin and mucous membranes. If reagents or specimens come in contact with sensitive area, wash with copious amounts of water.

Quality Control

Refer to BioGenex detection system documents for guidance on general quality control procedures.

Troubleshooting

Refer to the troubleshooting section in the documentation for BioGenex Detection Systems (or equivalent detection systems) for remedial actions on detection system related issues, or contact BioGenex Technical Support Department at 1-800-421-4149 or support@biogenex.com or your local distributor to report unusual staining.

Expected Results

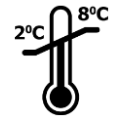





This antibody stains nucleus in positive cells in formalin-fixed, paraffin embedded tissue sections. An example image of a tissue section stained with this antibody can be found on the product page on the BioGenex website. Interpretation of the staining result is solely the responsibility of the user. Experimental results should be confirmed by a medically-established diagnostic product or procedure.

Limitations of the Procedure

Improper tissue handling and processing prior to immunostaining can lead to inconsistent results. Variations in embedding and fixation or the nature of the tissue may lead to variations in results. Endogenous peroxidase activity or pseudo peroxidase activity in erythrocytes and tissue biotin may result in non-specific staining based on the detection system employed. Tissues containing Hepatitis B Surface Antigen (HBsAg) may give false positive with horseradish peroxidase systems. Improper counterstaining and mounting may compromise the interpretation of results.

Bibliography

- Ruebner B H, et al. Cytomegalovirus infection. Viral ultrastructure with particular reference to the relationship of lysosomes to cytoplasmic inclusions. Am J Pathol 48:971-989, 1966.
- Niedobitek G, Finn T, Herbst H, Gerdes J, Grillner L, Landqvist M, Wirgart BZ, Stein H: Detection of cytomegalovirus by in situ hybridisation and immunohistochemistry using new monoclonal antibody CCH2: a comparison of methods. J Clin Pathol. 41:1005- 9 (1988).
- Ashley R, Peterson E, Abbo H, Gold D, Corey L: Comparison of monoclonal antibodies for rapid detection of cytomegalovirus in spin-amplified plate cultures. J Clin Microbiol. 27:2858-60 (1989).
- Zweyberg Wirgart B, Landqvist M, Hokeberg I, Eriksson BM, Olding-Stenkvist E, Grillner L: Early detection of cytomegalovirus in cell culture by a new monoclonal antibody, CCH2. J Virol Methods. 27:211-9 (1990).
- Canterino JE, McCormack M, Gurung A, Passarelli J, Landry ML, Golden M. Cytomegalovirus appendicitis in an immunocompetent host. J Clin Virol. 2016 May;78:9-11.

	Temperature Limitation		Batch Code
	Use By Date		Consult Instructions for Use
	Non-Sterile		Manufacturer

© 2020, BioGenex Laboratories. All rights reserved.

Category	Antibodies	Revision No.	D
Document No.	932-997M-EN	Release Date	22-Nov-2021